

eylys
Eye Care Kit

For Everyday Eye Comfort & Relaxation


A Cornea Specialist (FLVPEI) Initiative



AVAILABLE THROUGH JK MEDICAL IN LATIN AMERICA

Inside the Box

Directions for Use

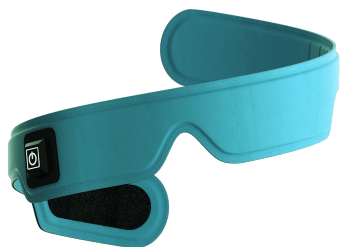

1 Warming eye device

2 Massaging eye roller

3 USB type-C to type-C cable

4 Welcome guide

DIRECTIONS FOR USE



Step 1: Switch on

Use the USB type-C to type-C cable to connect the Eylys warming eye device to a 5 volt USB type-C power source such as a mobile phone, power bank, or laptop.

Press the switch to turn on the eye device. It starts in Blue (less warm) mode. If needed, press the same button to change to Green (warm) or Red (more warm) mode, based on individual preference.



Step 2: Wear

Fasten the eye-contoured side of the device gently over closed eyes. Do not overtighten. Rest your head for better results. Once on, the device automatically switches off after 5 minutes.



Step 3: Massage

Use the massaging eye roller on closed eyes. Gently roll from top to bottom on the upper eyelid and from bottom to top on the lower eyelid. Repeat 5 times.

Features

Perfect fit. Portable comfort.

FEATURES

1 Perfect fit. Perfect pamper.

Eylys' patented 3D shape snugly fits your eyes to activate all upper and lower eyelid Meibomian glands.

2 Indulge. Anytime, anywhere.

Plug and power on the go with any USB type-C device, such as your phone, using a standard type-C to type-C cable.

3 3 temperature modes

With 3 safe and effective modes, patients can choose the one they find most soothing.

4 Quick and easy

Eylys' 5-minute therapy with automatic switch-off helps refresh your eyes in the blink of an eye.

5 Multi-level safety

Eylys includes multi-level safety and built-in thermal control systems.

6 Soft. Smooth. Soothing.

Made from hypoallergenic fabric with an antimicrobial coating, Eylys is super skin-friendly.

BENEFITS

- ☐ Eye Strain Relief
- ☐ Soothes Dry Eyes
- ☐ Deep Relaxation
- ☐ Natural Hydration
- ☐ Digital Protection

POSITIONING

- ☐ Medical
- ☐ Optical
- ☐ Pharmacy
- ☐ Corporate Gifting
- ☐ Hotels & Spas

Device Comparison: Eylys vs. Others

Eylys Eye Care Kit

	EYLYS	OTHERS
Design & USP	Designed by a Cornea Specialist (FLVPEI). 3D eye-fit technology* matches every unique eye contour for best-in-class comfort and relief.	Bulky granules / gel mask (if non-electric) or flat mask (if electric) do not touch the eyelid skin completely, leading to incomplete relief.
Features	Massaging eye roller. Powered by USB-C phone. 3 temperature modes: all safe and effective. Auto-off in just 5 minutes.	If non-electric: none. If electric: very high (unsafe) or low (not effective) temperature modes, with auto-off at up to 60 minutes.
Additional Requirements To Use Device	Uses a standard USB-C mobile cable that can be plugged into a USB-C phone to power the device.	Needs a microwave, gas stove, or electric kettle (if non-electric mask), or a proprietary cable (if electric mask).
Thermal Control & Safety	Consistent performance despite changes in ambient temperature. Uniform heat distribution over the entire upper and lower eyelids, with no hot or cold spots.	If non-electric: very high initial temperature (likely to cause thermal injury), followed by a fast drop to ineffective levels. If electric: slow temperature rise with non-uniform heat distribution, causing hot and cold spots.
Backup safety if device is abused	Auto-detection of internal circuit damage and permanent switch-off.	None.

*Patent-pending. All intellectual properties are owned by or licensed to Oculab Pvt. Ltd., India.



EYE CARE KIT

Consistent warmth. Simple routine.



AVAILABLE THROUGH JK MEDICAL IN LATIN AMERICA

Availability and authorized product claims vary by country.

JK Medical, Inc.

Latin America Distributor Contact

Office: +1 (352) 324-3091

Mobile: +1 (407) 342-8575

Email: info@jkmedicalinc.com

Web: www.jkmedicalinc.com

Manufacturer

Oculab Private Limited

Mumbai, India

www.oculab.in

For personal eye-care and wellness use. This brochure is a commercial overview and does not replace the user manual or professional medical advice. Product availability, regulatory status, labeling, indications, warnings, packaging, and included accessories vary by country. Confirm the locally authorized version before sale or use.